

# Blood Purification Products




## Standard Kits include:

1. Hemodialysis Catheter
2. Vessel Dilator
3. Introducer Needle
4. Syringe
5. Guide-wire
6. Adhesive Wound Dressings
7. Heparin Caps
8. Scalpel
9. Needle with Suture

## Optional Compound Kits include:

1. Hemodialysis Catheter Standard Kit Accessories
2. 5ml Syringe
3. Surgical Gloves
4. Surgical Pledget
5. Surgery Sheet
6. Surgery Towel
7. Sterile Brush
8. Gauze Pad

# Hemodialysis Catheter

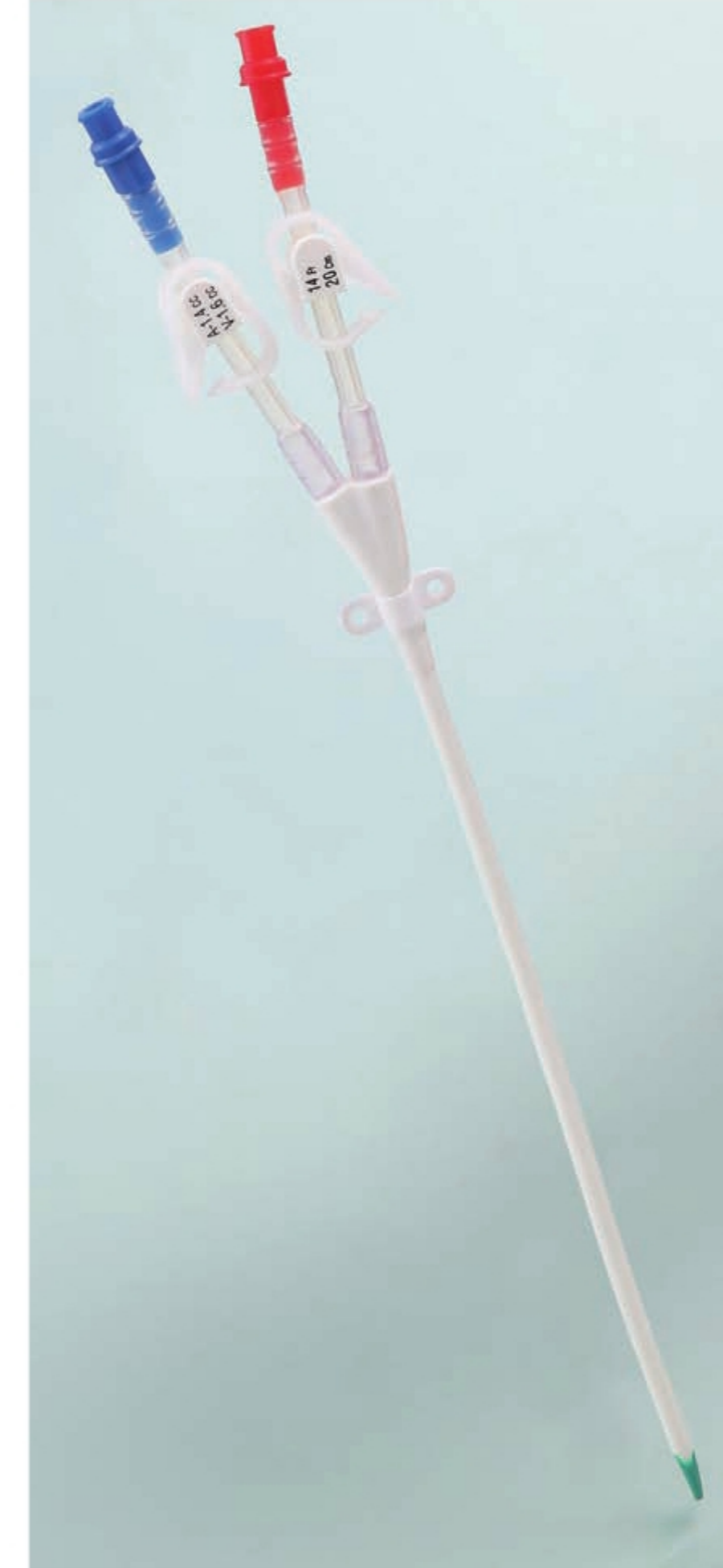

## Soft Tip

Tapered tip facilitates ease of insertion without use of peel-away sheath. Reduces vessel trauma during insertion.



## Side Holes:

2.5cm far from end hole reduces recirculation and wall-suction. Minimizes risk of clot formation.



## Radiopaque

Radiopaque, facilitates quick visualization under X-Ray.



## Rotatable Suture Wing

Facilitates skin inspection. Avoid blood vessel wall-suction. Minimizes risk of exit infections. No tube-kinking occurred.



## Silicone Extension Tube

Increases patient comfort, will not crimp over time, able to visualize fluids.



## Multi-lumen available

Single, double, triple lumen.



## Compound-packaging Tray

Various Accessories for choice.



## Made from PU (medical grade)

Better flexibility and anti-kinking performance.

## Pre-curved Type

More comfortable on patients.





## Order Information

Product Code	Lumen	Catheter Size	Catheter Length(cm)	Dilator	Guidewire	Shape of Extension Tube
FH-1721	Single	7Fr	5	8F×10cm	0.035"×50cm	Straight
FH-1721W	Lumen					Curved
FH-1722	Single	7Fr	8	8F×10cm	0.035"×50cm	Straight
FH-1722W	Lumen					Curved
FH-1723	Single	7Fr	10	8F×10cm	0.035"×50cm	Straight
FH-1723W	Lumen					Curved
FH-1724	Single	7Fr	13	8F×10cm	0.035"×50cm	Straight
FH-1724W	Lumen					Curved
FH-1725	Single	7Fr	16	8F×10cm	0.035"×50cm	Straight
FH-1725W	Lumen					Curved
FH-1726	Single	7Fr	20	8F×10cm	0.035"×60cm	Straight
FH-1726W	Lumen					Curved

Product Code	Lumen	Catheter Size	Catheter Length(cm)	Dilator	Guidewire	Shape of Extension Tube
FH-1821	Single	8Fr	5	9F×10cm	0.038"×50cm	Straight
FH-1821W	Lumen					Curved
FH-1822	Single	8Fr	8	9F×10cm	0.038"×50cm	Straight
FH-1822W	Lumen					Curved
FH-1823	Single	8Fr	10	9F×10cm	0.038"×50cm	Straight
FH-1823W	Lumen					Curved
FH-1824	Single	8Fr	13	9F×10cm	0.038"×50cm	Straight
FH-1824W	Lumen					Curved
FH-1825	Single	8Fr	16	9F×10cm	0.038"×50cm	Straight
FH-1825W	Lumen					Curved
FH-1826	Single	8Fr	20	9F×10cm	0.038"×60cm	Straight
FH-1826W	Lumen					Curved

## Order Information

Product Code	Lumen	Catheter Size	Catheter Length(cm)	Dilator	Guidewire	Shape of Extension Tube
FR/FH-2621	Double	6.5Fr	5	7F×10cm	0.021"×50cm	Straight
FR/FH-2621W	Lumen					Curved
FR/FH-2622	Double	6.5Fr	8	7F×10cm	0.021"×50cm	Straight
FR/FH-2622W	Lumen					Curved
FR/FH-2623	Double	6.5Fr	10	7F×10cm	0.021"×50cm	Straight
FR/FH-2623W	Lumen					Curved
FR/FH-2624	Double	6.5Fr	13	7F×10cm	0.021"×50cm	Straight
FR/FH-2624W	Lumen					Curved
FR/FH-2625	Double	6.5Fr	16	7F×10cm	0.021"×50cm	Straight
FR/FH-2625W	Lumen					Curved
FR/FH-2626	Double	6.5Fr	20	7F×10cm	0.021"×50cm	Straight
FR/FH-2626W	Lumen					Curved
FR/FH-2921	Double	8.5Fr	5	9F×10cm	0.025"×50cm	Straight
FR/FH-2921W	Lumen					Curved
FR/FH-2922	Double	8.5Fr	8	9F×10cm	0.025"×50cm	Straight
FR/FH-2922W	Lumen					Curved
FR/FH-2923	Double	8.5Fr	10	9F×10cm	0.025"×50cm	Straight
FR/FH-2923W	Lumen					Curved
FR/FH-2924	Double	8.5Fr	13	9F×10cm	0.025"×50cm	Straight
FR/FH-2924W	Lumen					Curved
FR/FH-2925	Double	8.5Fr	16	9F×10cm	0.025"×50cm	Straight
FR/FH-2925W	Lumen					Curved
FR/FH-2926	Double	8.5Fr	20	9F×10cm	0.025"×50cm	Straight
FR/FH-2926W	Lumen					Curved
FR/FH-2124	Double	11.5Fr	13	11.5F×15cm	0.038"×50cm	Straight
FR/FH-2124W	Lumen					Curved
FR/FH-2125	Double	11.5Fr	16	11.5F×15cm	0.038"×50cm	Straight
FR/FH-2125W	Lumen					Curved
FR/FH-2126	Double	11.5Fr	20	11.5F×15cm	0.038"×60cm	Straight
FR/FH-2126W	Lumen					Curved
FR/FH-2127	Double	11.5Fr	25	11.5F×15cm	0.038"×60cm	Straight
FR/FH-2127W	Lumen					Curved

Product Code	Lumen	Catheter Size	Catheter Length(cm)	Dilator	Guidewire	Shape of Extension Tube
FR/FH-2224	Double	12Fr	13	12F×15cm	0.038"×50cm	Straight
FR/FH-2224W	Lumen					Curved
FR/FH-2225	Double	12Fr	16	12F×15cm	0.038"×50cm	Straight
FR/FH-2225W	Lumen					Curved
FR/FH-2226	Double	12Fr	20	12F×15cm	0.038"×60cm	Straight
FR/FH-2226W	Lumen					Curved
FR/FH-2524	Double	14Fr	13	14F×15cm	0.038"×50cm	Straight
FR/FH-2524W	Lumen					Curved
FR/FH-2525	Double	14Fr	16	14F×15cm	0.038"×60cm	Straight
FR/FH-2525W	Lumen					Curved
FR/FH-2526	Double	14Fr	20	14F×15cm	0.038"×70cm	Straight
FR/FH-2526W	Lumen					Curved
FR/FH-3224	Triple	12Fr	13	12F×15cm	0.038"×50cm	Straight
FR/FH-3224W	Lumen					Curved
FR/FH-3225	Triple	12Fr	16	12F×15cm	0.038"×50cm	Straight
FR/FH-3225W	Lumen					Curved
FR/FH-3226	Triple	12Fr	20	12F×15cm	0.038"×60cm	Straight
FR/FH-3226W	Lumen					Curved

### Remarks:

1. The product codes in this chart are for HC kits with straight introducer needle configuration.
2. HC kits with Y-shaped introducer needle configuration are also available upon discussion. The product code of HC kit with Y-shaped introducer needle configuration is different from the product code of HC kit with straight introducer needle configuration. Before you order, please clarify which configuration you request, either straight or Y-shaped introducer needle.
3. Pediatric types of HC include double lumen 6.5Fr and 8.5Fr.
4. Adult types of HC include single lumen 7Fr, 8Fr ; double lumen 11.5Fr, 12Fr, 14Fr and triple lumen 12Fr.
5. "FR" in the product code means with soft tip.
6. "FH" in the product code means without soft tip (i. e. with relatively harder tip).
7. Pre-curved catheter is available for adult type in double lumen 11.5Fr, 12Fr, 14Fr.