



DISCOVER
HEALTHCARE
EXPERTISE



HEALTHCARE DOSSIER

Godrej interio®



HEALTHCARE AND US

Health and care are inseparable. Health improves with the right care. Health related caregivers have, thus, the noblest profession in the world. It is to support and enhance the work of these tireless medics that Godrej Interio has developed solutions to all the problem areas within this arena.

SHAPED BY CHALLENGES

In an environment where every day is a challenge, our products act as facilitators of care and as a support system that supplements treatments. They deliver significant benefits to the patient and caregiver alike. It is these very challenges that medics face on a day-to-day basis that have been instrumental in providing inputs leading to the design and development of this intuitive, dynamic and ergonomic healthcare product range.



DESIGNED WITH CARE, EXPERIENCE AND EXPERTISE

Further to simply improving upon established solutions trusted globally by the medical fraternity, our design teams concentrate on creating products that cater to key customer needs common to a vast variety of healthcare environments – so that the products they create have multiple-use, yet are specialized for specific tasks.

These products are made keeping in mind the following:

Unparalleled Safety as per Global Standards and Certifications:

In line with strict International norms, our products are created to have no sharp edges, no electrostatic charge, no cold shock, all entrapment zones taken care of, fall protection, adherence to the most stringent global standards and quality certifications.



Comfort and Care for both patients and caregivers:

Easy-to-use controls/crank handles, ergonomic design, easy maneuverability and compatibility with associated products and surroundings are some of the features which make certain that there is no compromise on comfort and care for all concerned.

Hygiene & Maintenance:

Our products carry a bacteriostatic powder coating that ensures effective infection control, better pre-treatment procedure and simple construction with easy access for cleaning and fumigation.



GODREJ INTERIO: HOLISTIC DESIGN AND BUILD SOLUTIONS

⌘ Creation of environments that work in harmony towards patient recovery, incorporating a complete healthcare system.



⌘ Strong emphasis on creating facilities that promote well-being through interior designs that emphasize professionalism and medical expertise as well as bring about an international feel through improved workflow patterns and processes.



⌘ Design team with years of practical experience and exposure to every Interior, HVAC, electrical, and plumbing design – with experience in working processes as builders and installers.



Godrej Interio offers you its entire basket of Healthcare as well as other application area furniture. Ease of use, safety and ease of maintenance are the underlying benefits of this range, allowing us to be a one-stop solution for all your needs.

MOTORIZED BEDS

CHRYSLIS

- 1

SIDEBOARDS:
Shelters the patient from falls
in all bed positions.
- 2

CAREGIVER
CONTROL PANEL:
Ensures caregiver efficiency
and patient safety, with
enhanced safety features.
- 3

ENCLOSED FORM OF
HEAD AND LEG BOARD:
Aligned with sideboard
for enhanced protection.
- 4

BUMPERS:
Absorbs impact during mobility.
- 5

LINEN HOLDER:
Facilitates changing of bed linen.
- 6

CASTOR BRAKES WITH
DIRECTIONAL LOCK AND
AUXILIARY BRAKING SYSTEM:
For guaranteed safety and
control during mobility.
- 7

DROP DOWN
SIDEBOARD MECHANISM:
Controlled, dampened release
mechanism for safe and silent
sideboard movements.



- PATIENT
CONTROL PANEL:
Gives autonomy to the patient.

8
- UNCLUTTERED UNDER
STRUCTURE:
For clean and sleek look and
ease of maintenance.

9
- WIRE WINDER:
For efficient wire management
and clutter-free access.

10
- ANGLE INDICATOR:
Gives clarity and precision in adjustment
of back-rest, knee-rest and Trendelenburg/
Reverse Trendelenburg positions.

11
- TWIN WHEEL CASTORS
WITH CENTRAL LOCKING:
Provides ease of mobility and stability.

12
- CPR LEVER:
Helps achieve flat
position of backrest.

13
- LEG ELEVATION
ON RATCHET:
Helps adjust leg rest section.

14
- REMOVABLE MATTRESS
PLATFORM:
ABS removable surface made easy
in cleaning and fumigation.

15
- REMOVABLE HEAD/FOOT BOARD
Head/Foot board is firmly fixed to the bed.
Easy removability for access to the patient.

16

WORLD CLASS STANDARDS

IEC 60601-1	Medical Electrical Equipment: General requirements for safety and essential performance.
IEC 60601-2-52*	Medical Electrical Equipment: Particular requirements for safety and essential performance of medical beds.
EN 1970:2000	Adjustable bed for disabled persons. (Requirements & Tests)
IEC 60601-1-2	Electromagnetic compatibility. (Requirements & Tests)
IEC 60601-1-4	General Requirements for Safety: Programmable electrical medical systems.
ISO 14971	Medical Devices: Application of risk management to medical devices.
EN 980:2008	Symbols for use in the labelling of medical devices.

*IV Pole available

Chrysalis is a range of 5 function, fully motorized patient care beds comprising of 24 variants, which is designed for care of patients in the spectrum of Intensive Care to Post Operative Care. It is suitable for ICCUs, ICUs and Private Rooms / Dual Occupancy Rooms.

CHRYsalis NEO

MOTORIZED BEDS

Chrysalis Neo is a range of fully motorized patient care beds comprising of 16 variants, which is designed for care of patients in the spectrum of Intensive Care to Post Operative Care. It is suitable for ICCUs, ICUs and Private Rooms / Dual Occupancy Rooms of Private and Corporate Hospitals.




CHRYsalis NEO ICU XP
with CPR

Castor Options	Castors
Twin Wheel Central Locking	
Single Wheel Central Locking	
Single Wheel Diagonal Locking	


Colour Options
Sea Green
Peach
Grey

Dimensions (In mm)			
Size	L	W	H
	2206	1010	450 to 825


ADDITIONAL ACCESSORIES WITH SELECTED CHRYsalis BEDS

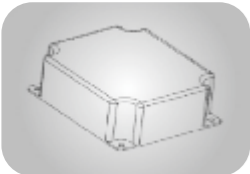


Battery Backup

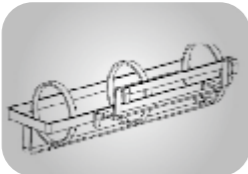

Cylinder Holder

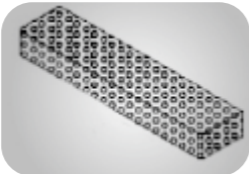

Mattress Extension

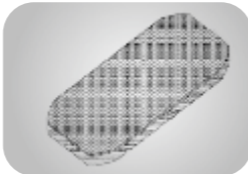

Chrysalis Mattress

CHRYsalis

MOTORIZED BEDS




CHRYsalis ICU XP
with CPR

Castor Options	Castors
Twin Wheel Central Locking	
Single Wheel Central Locking	
Single Wheel Diagonal Locking	

Colour Options
Sea Green
Peach
Grey

Dimensions (In mm)			
Size	L	W	H
	2206	1010	450 to 825

* All accessories shown are optional and available at extra cost

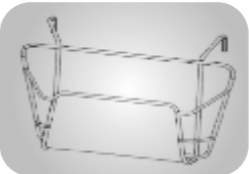

Document Holder



X-RAY CASSETTE



Lite IV Pole



SEMI-MOTORIZED BEDS

E-VITA

1 REMOVABLE HEAD/FOOT BOARD
Head/Foot board are firmly fixed to the bed.
Easily removable for access to the patient.


2 CAREGIVER CONTROL PANEL:
Ensures caregiver efficiency and patient safety, with enhanced safety features.

3 SWIVEL SIDEBOARD MECHANISM:
Easy to operate side boards with a Self Locking System (SLS) at 180°.

4 BUMPERS:
Absorbs impact during mobility.

5 UNCLUTTERED UNDER STRUCTURE:
For clean and sleek look and ease of maintenance.

6 CASTOR BRAKES WITH DIAGONAL LOCK
For guaranteed safety and control during mobility.

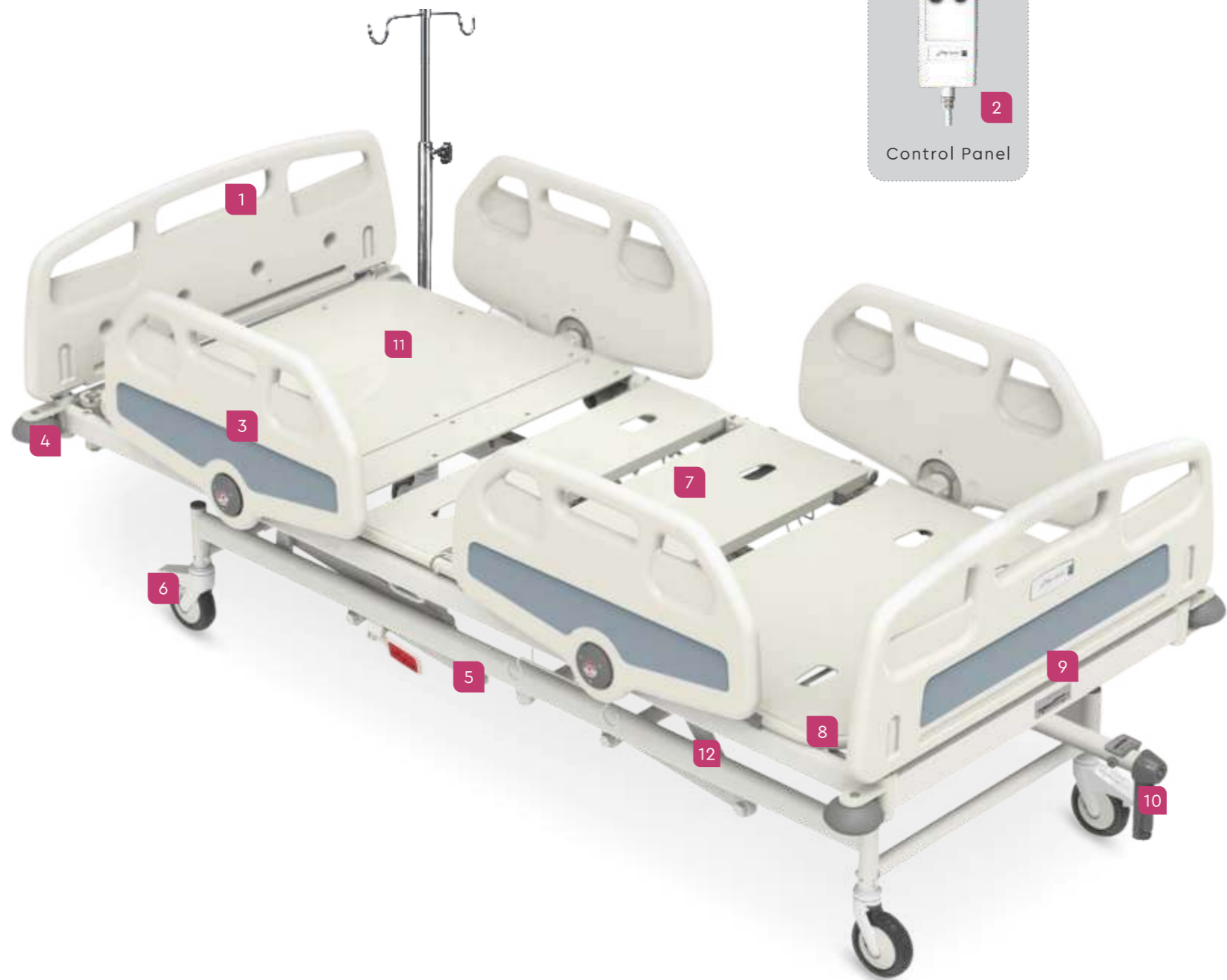

Control Panel

7 METAL LYING SURFACE
Metal lying surface with mattress retainer.

8 LOWER LEG REST ELEVATION ON RATCHET:
Helps adjust leg rest section.

9 GRAPHICAL SYMBOLS
Indicate position adjustments achieved by corresponding handles, for clarity and simplicity of use.

10 CRANK HANDLE
Ergonomically designed. Features outward locking in the 90° position. Minimizes pressure on the wrist.

11 X-RAY MECHANISM
A top-loading SS plate for loading X-ray cassette, locked in place by spring loaded knobs.


12 SLEEK BASE FRAME
Sleek frame with large clear space and high ground clearance with high degree of stability.

E-Vita is a range of 2 or 3 function motorized Patient Care Beds, with 8 variants. It is suitable for ICUs and Private Rooms/Dual Occupancy Rooms of Private, Corporate and Government Hospitals.




SEMI-MOTORIZED BEDS

E-VITA ICU XP Bed
With Manual Trendelenburg

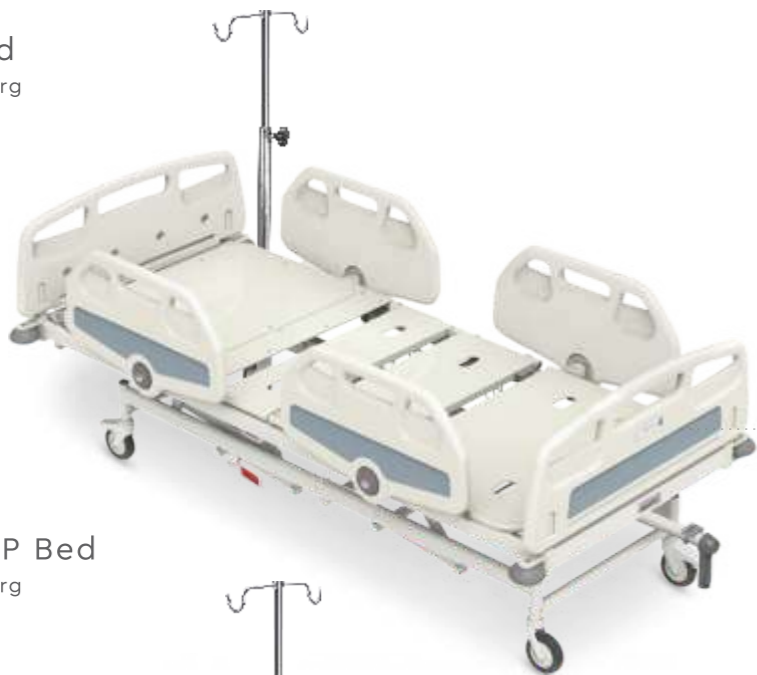

E-VITA LITE ICU XP Bed
With Manual Trendelenburg

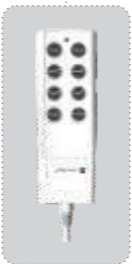

- Colour Options
- Sea Green
 - Peach
 - Grey

Dimensions (In mm)

Size	L	W	H
	2312	1060	490 to 770

Motorized	●	●			●
Manual			●	●	

E-VITA ICU XP Bed
With Manual Trendelenburg
and Knee rest




E-VITA LITE ICU XP Bed
With Manual Trendelenburg
and Knee rest



Motorized	●				●
Manual		●	●	●	


ADDITIONAL ACCESSORIES WITH SELECTED CHRYSALIS BEDS

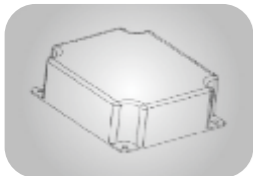

Battery Backup




Lite IV Pole

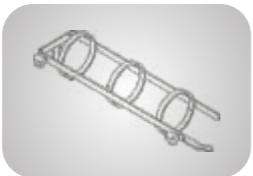

Cylinder Holder



Document Holder

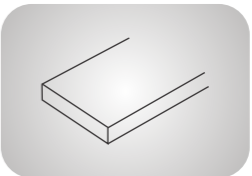

Suitable Mattress



X-Ray Cassette

* All accessories shown are optional and available at extra cost



MANUAL BEDS

VITA

1 REMOVABLE HEAD/FOOT BOARD:
Head/Foot board are firmly fixed to the bed.
Easily removable for access to the patient.

2 SLEEK BASE FRAME:
Sleek frame with large clear space and high ground clearance with high degree of stability.

3 SWIVEL SIDEBOARD MECHANISM:
Easy to operate sideboards with a Self Locking System (SLS) at 180°.

4 BUMPERS:
Absorbs impact during mobility.

5 UNCLUTTERED UNDER STRUCTURE:
For clean and sleek look and ease of maintenance.

6 CASTOR BRAKES WITH DIAGONAL LOCK:
For guaranteed safety and control during mobility.

7 HEIGHT INDICATOR:
Trendelenburg operating height indicator marks the correct bed height for adjusting Trendelenburg position.

8 LOWER LEG REST ELEVATION ON RATCHET:
Helps adjust leg rest section.

9 GRAPHICAL SYMBOLS:
Indicate position adjustments achieved by corresponding handles, for clarity and simplicity of use.

10 CRANK HANDLE:
Ergonomically designed. Features outward locking in the 90° position. Minimizes the pressure on the wrist.

11 X-RAY MECHANISM:
A top-loading SS plate for loading X-ray cassette, locked in place by spring loaded knobs.

12 METAL LYING SURFACE:
Metal lying surface with mattress retainer.



Vita is a range of manual ICU beds for use in Pre and Post-operative Care Rooms of Government and Tertiary Care Hospitals.



VITA

MANUAL BEDS



VITA ICU Bed
With Metal Surface



Colour Options

- Sea Green
- Peach
- Grey

Dimensions (In mm)


Size	L	W	H
	2207	1060	485 to 750



ADDITIONAL ACCESSORIES WITH SELECTED CHRYSALIS BEDS



Lite I.V. Pole

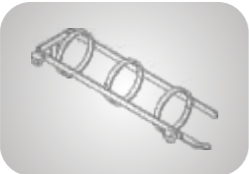

Cylinder Holder



X-ray Cassette



VITA LITE

MANUAL BEDS



VITA Lite ICU



Colour Options

- Sea Green
- Peach
- Grey

Dimensions (In mm)


Size	L	W	H
	2207	1060	485 to 750



* All accessories shown are optional and available at extra cost



Crank Handle

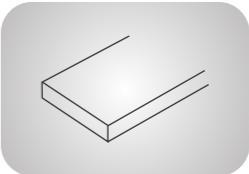

Suitable Mattress



Document Holder



AYUSH PLUS

MANUAL BEDS



AYUSH PLUS Fowler*

Ayush Plus Fowler Bed is a Mechanical Two Function Bed which is designed for Patient Care in non-critical settings. These beds are perfect especially for hospitals dealing with high occupancy, rough usage and seeking a simple robust product.



Size	L	W	H
	2143	915	560





AYUSH PLUS Semi-Fowler*



Size	L	W	H
	2143	915	560






Ayush Plus Sheet Bed*

Size	L	W	H
	1998	915	560



*SS PLT Removable Head/Foot-board

ADDITIONAL ACCESSORIES WITH SELECTED CHRYSALIS BEDS

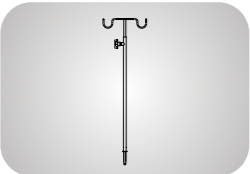

Lite I.V. Pole




Crank Handle



Side Rail



AYUSH PLUS Fowler PP#



Size	L	W	H
	2141	990	560





AYUSH PLUS Semi-Fowler PP#



Size	L	W	H
	2141	990	560



#PP Removable Head/Foot-board



Ayush Plus Sheet PP#


Size	L	W	H
	1985	990	565



* All accessories shown are optional and available at extra cost



Document Holder

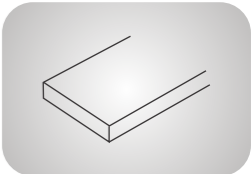


Suitable Mattress



AYUSH

MANUAL BEDS


AYUSH Fowler



Size	L	W	H
	2222	922	610




AYUSH Semi-Fowler



Size	L	W	H
	2222	922	610




AYUSH Sheet Bed

Size	L	W	H
	2030	915	614




AYUSH Attendant Bed

Size	L	W	H
	2025	774	550



Size	L	W	H
	1980	700	100




*Mattress of suitable size available for Attendant Bed

ADDITIONAL ACCESSORIES WITH SELECTED CHRYSALIS BEDS



Lite I.V. Pole

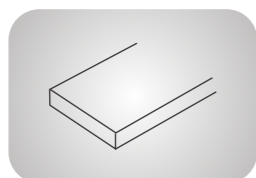

Suitable Mattress



Crank Handle



Side Rail



Document Holder

* All accessories shown are optional and available at extra cost

EXAMINATION TABLE

ELIXIR

1 BACKREST ON GAS-SPRING:
Headrest can be adjusted with a gentle press to the U-frame handle. The easy to adjust headrest can be adjusted up to 30 degrees.

2 TISSUE ROLL HOLDER:
Tissue roll holder enables fresh sheets of disposable tissues providing each patient a fresh hygienic surface.

3 BP TRAY:
Compact tray tucked below the mattress surface and can be swivelled out when needed.

4 TRAPEZOID DESIGN:
Narrow base with high stability and ample comfort. Offers clear standing room and optimum patient access to the caregiver.

5 INTEGRATED UPHOLSTERY WITH ROUNDED CORNERS:
Soft rounded corners ensure patient safety and provide much needed comfort.

6 LEVELLING SCREWS:
Enhances product stability over uneven surfaces.

7 FOOTSTEP:
Footstep with anti-skid rubber mat ensures safe patient mount-dismount.



 ELIXIR EXAMINATION COUCH
FULL STORAGE

The ELIXIR range of examination couches and tables designed for OPD consultation rooms have an in-built array of innovative, value-added features that ensure superior patient comfort, hygiene and ease of use.



ELIXIR

EXAMINATION TABLES



ELIXIR Examination Couch
Compact Storage



ELIXIR Examination Table
Drawer Storage



ELIXIR Examination Table
Basic



Colour Options for Mattress Upholstery

Beige Maroon Dark Grey

Dimensions (In mm)

Size	L	W	H
	1957	625	810



ELIXIR

GYNAEC TABLES

The ELIXIR range of Gynaecological Examination Couches and Tables is designed for Maternity Hospitals, Nursing Homes and Gynaec/Maternity/Obstetric Wings of large Private and Corporate Hospitals wherein Gynaecological consultations are carried out. It is uniquely built for a superior degree of patient comfort and hygiene level.



ELIXIR Gynaec
Examination Couch
With Drawer Storage



ELIXIR Gynaec
Examination Couch
Without Drawer Storage

Colour Options for Mattress Upholstery

Beige Maroon Dark Grey

Dimensions (In mm)

Size	L	W	H
	1954	710	810

EXAMINATION TABLES

Examination Table
With Ratchet Backrest


Size	L	W	H
	1830	600	850



Examination Table
Gynaec

Size	L	W	H
	2060	850	850

OPTIONAL ACCESSORIES

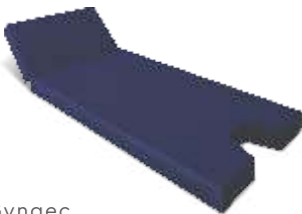

Mattress for Gynaec
Examination Table



Mattress for
Examination Table
on Ratchet

* All accessories shown are optional and available at extra cost

PAEDIATRICS & INFANT CARE

Infants Crib

Size	L	W	H
	1008	536	895



Paediatric Bed

Size	L	W	H (till Bed Surface)
	1572	839	1218



Crib Arrangement on Bed

Size	L	W	H
	908	416	784




Mattress for Paediatric Bed

Size	L	W	H
	1415	785	75



Mattress for Crib

Size	L	W	H
	615	315	50

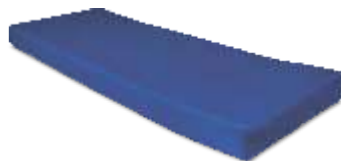

*Mattress available at extra cost

APTIS / AYUSH

LOCKERS

APTIS Locker

With Double Drawer and Cabinet

Size	L	W	H
	490	410	880



APTIS Locker

With Drawer and Cabinet

Size	L	W	H
	490	410	941



APTIS Locker

With Cabinet

Size	L	W	H
	490	410	941



Colour Options: Sea Green Peach Grey

AYUSH Delux Locker

With Drawer and Cabinet

Size	L	W	H
	421	402	840



AYUSH Simple Locker

With Cabinet

Size	L	W	H
	404	402	835



APTIS / AYUSH

OVERBED TABLES

ADEO Overbed Table

Height Adjustable With Gas Lift

Size	L	W	H
	896	395	792 to 1092



APTIS Overbed Table

On Gaslift

Size	L	W	H
	896	395	792 to 1092



APTIS Overbed Table

On Crank Lift

Size	L	W	H
	899	555	970 to 1170



AYUSH Overbed Table

Fixed Height

Size	L	W	H
	1050	380	950



PATIENT MOBILITY PRODUCTS_____

 Rapid Neo Emergency Recovery Trolley



Size	L	W	H
	2139	936	700 to 1080



*MATTRESS OPTION AVAILABLE AT COST

 Trauma Care Trolley
Height Adjustable



Size	L	W	H
	1900	750	520 to 750



 Stretcher on Trolley

Size	L	W	H
	2005	666	830

MATTRESS OPTION AVAILABLE AT COST.




 Folding Wheelchair With Cushion

Size	L	W	H
	790	600	870



TROLLEYS & CARTS

 Instrument Trolley
With Bowl and Bucket

Size	L	W	H
	1156	531	915



 Instrument Trolley

Size	L	W	H
	775	531	915



 Crash Cart

Size	L	W	H
	1048	475	1555



 Medicine Trolley

Size	L	W	H
	902	532	915



 Mayo Trolley

Size	L	W	H
	504	650	814 to 1340



 Instrument Cabinet

Size	L	W	H
	900	450	1800



 Oxygen Cylinder Trolley

Size	L	W	H
	410	388	1112



SS & MS ACCESSORIES

Basin Stand

Size	L	W	H
	425	425	850



Kick Bucket

Size	L	W	H
	517	517	210



Square Stool

SS, Top

Size	L	W	H
	411	411	522



Round Stool

SS, Top

Size	L	W	H
	538	538	470 to 655



Single Footstep

Size	L	W	H
	485	335	210



Double Footstep

Size	L	W	H
	665	485	393



Three Fold Partition Screen

Size	L	W	H
	2637	640	1720



Soiled Linen Trolley

Size	L	W	H
	625	625	955



SS & MS ACCESSORIES

Saline Stand
SS, 5 Castors

Size	W (base diameter)	H
	652	1655 to 2365



Saline Stand
MS, 4 Castors

Size	W (base diameter)	H
	695	2121



Infusion Pole

Size	W (base diameter)	H
	366	1007 to 1300 From Bed Surface Top



Lite I.V. Pole

Size	W (base diameter)	H
	250	743 to 1150 From Bed Surface Top



Railings
MS & SS



Size	L	W	H (from Bed surface)
	1473	89	516

MATTRESSES

CHRYSLIS

Suitable for all Chrysalis
and Chrysalis Lite Beds

Size	L	W	H
	1905	840	100



4 Section

Suitable for E-vita, E-vita Lite, Vita,
Vita Lite ICU Beds, Vita Lite Fowlers,
Ayush Plus Fowlers, Ayush Fowlers

Size	L	W	H
	1965	830	100



2-Section

Suitable for Vita Lite,
Ayush Plus, Ayush Semi-Fowlers


Size	L	W	H
	1965	830	100



Single Section

Suitable for Ayush Plus,
Ayush Sheet Beds

Size	L	W	H
	1900	900	75,100

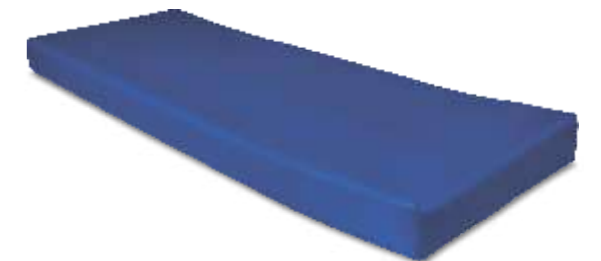

*CRIB/EXAMINATION/ATTENDANT/ERT



The entire gamut of services that can be offered under the Godrej Interio umbrella will include:



MEDICAL HEADS / SENIOR DOCTORS ROOM



Impress



Numero Uno



Maestro



Lead



NRG Plus



Pulse



TRAINING ROOM



Fold



Unlearn



Scintilla
With Desklet



Namaste



RECORD ROOM

Optimizer



Personal Locker Unit

Gain



LIBRARY

Book Rack

Glass Door Storewell

Storewell Plain



MEETING ROOM

Conclave



ADMIN AND NURSING STATION

Wish



Spacio



LFC



VSDU



CAFETERIA

Time Out



Cantina



Café Pod



MPC - 4



LAB

Aeolus



Fusion





OPDs



Gallop



Enterprise



Work



WAITING AREA



Matrix



Pisa Table



IMPECCABLE CREDENTIALS




CII DESIGN EXCELLENCE AWARD



INDIA DESIGN MARK



USFDA CERTIFICATE

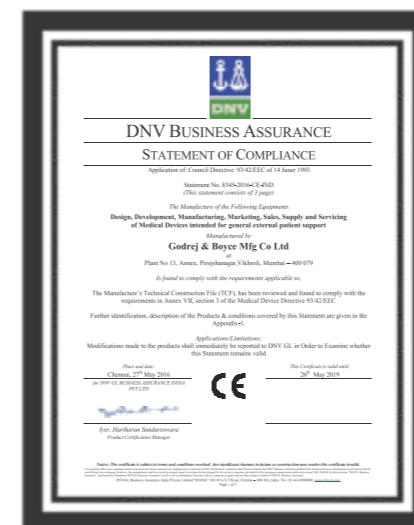

CE CERTIFICATE



ISO 13485



INDIA DESIGN MARK



AVAILABLE IN OVER 600 LOCATIONS ACROSS INDIA

Put some play into your work by visiting



An initiative by Godrej Interio, the first Indian furniture company with **GREENGUARD** certified products

MODULAR FURNITURE • DESKING • SEATING • STORAGE • INTERIORS • LAB • MARINE • HEALTHCARE • HOME

Local Address

National

North: Delhi, Chandigarh, NCR, Lucknow

South: Chennai, Bengaluru, Hyderabad, Kochi

West: Mumbai, Ahmedabad, Bhopal, Pune

East: Kolkata, Bhubaneswar, Guwahati, Patna

In keeping with the Godrej policy of continuous development and improvement, the design, dimensions and specifications may be updated without prior notice or obligations. Form, Colour and Finish of product in this document are representative of design intent and may vary from actual. Colour options depicted in this catalogue are as close to what modern printing processes allow. All dimensions are in 'cm' and are subject to commercial tolerance. Our products are meant for indoor use only.

Follow us on



/godrejinterio



/godrej-interio



/godrejinterio



/godrejinterio4u